

Eliminate Plate Bending

CNC Milled Titanium

Customized Plate Design

Optimized Surgical Outcomes

Surgeons choose **AccuPlate®** to optimize surgical outcomes, save time in the OR, and remove mechanical stress induced by manual plate bending.


Indications for Use

The MedCAD® AccuPlate® Patient-Specific Plate is intended for prescription use in oral and maxillofacial surgery, trauma and reconstructive surgery.


Specific Indications for Use

- Primary mandibular reconstruction with bone graft
- Temporary bridging until delayed secondary reconstruction
- Secondary mandibular reconstruction
- Comminuted mandibular fractures
- Fractures of edentulous and/or atrophic mandibles
- Unstable mandibular fractures
- Maxillary reconstruction with or without bone graft
- Maxillary trauma


Condyle-to-Condyle

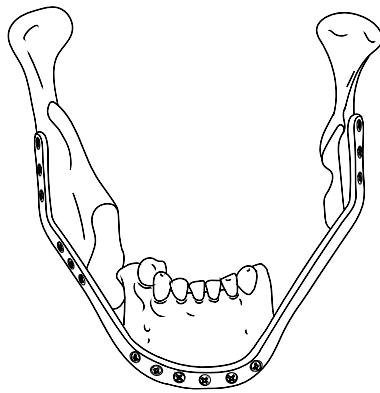


Angle-to-Angle

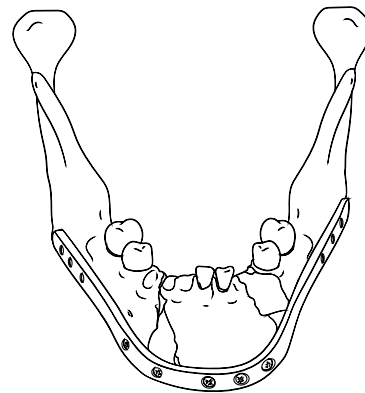

Hemi Right

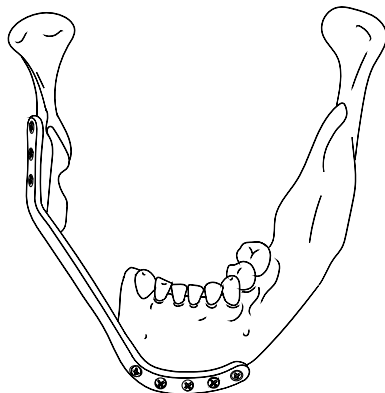


Hemi Left

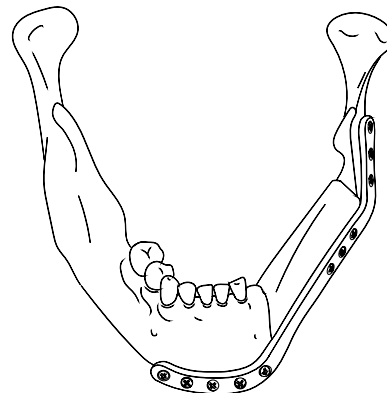

FEATURES

ACCUPLATE® Patient-Specific Plate

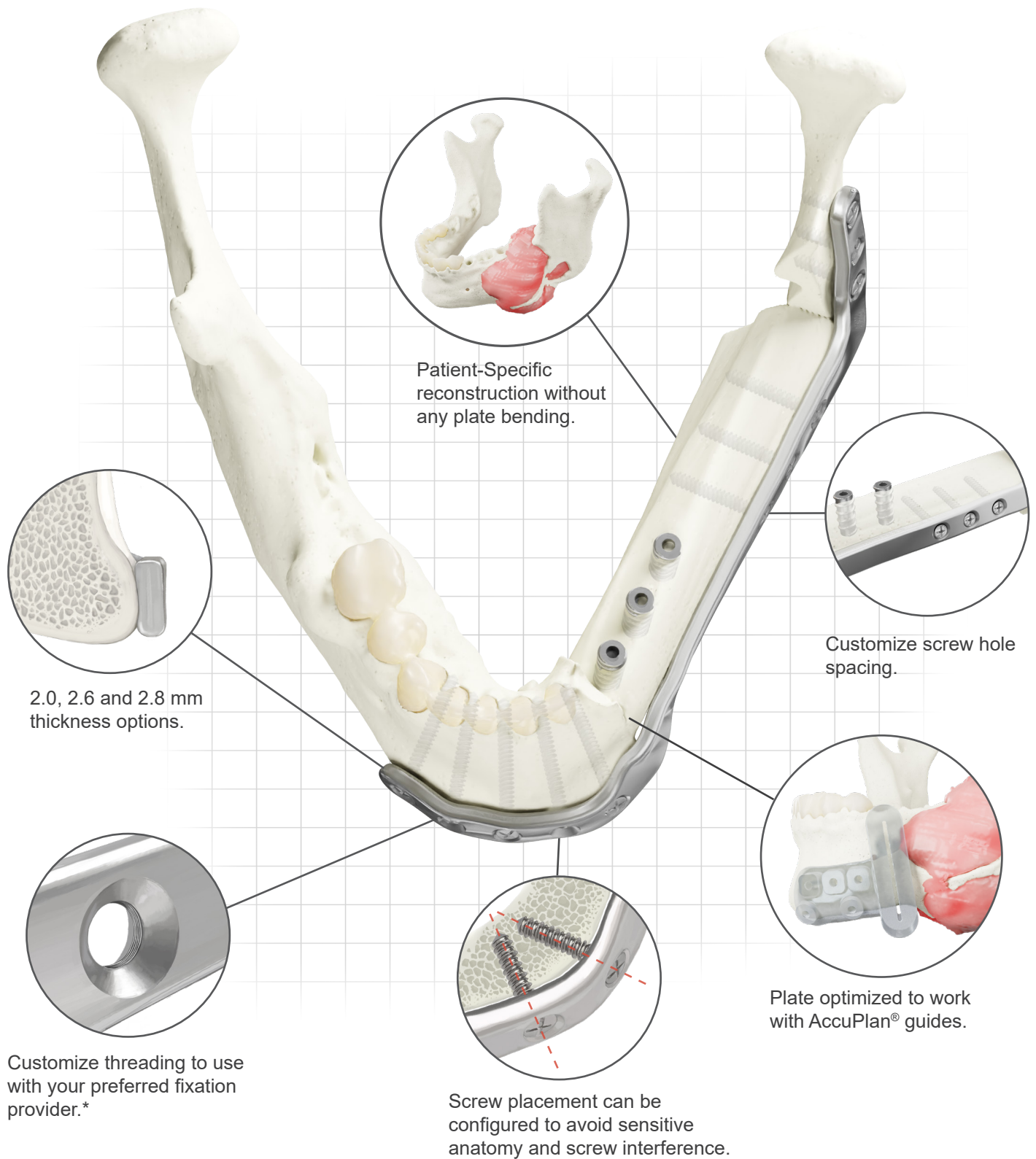

**Compatible with MedCAD-approved fixation options.*

Web
MedCAD.com/Start

Phone
+1 (214) 453-8864

MedCAD
501 S 2nd Ave, Ste A-1000
Dallas, TX 75226



ACCUPLATE® Patient-Specific Plate

Benefits

Patient-Specific plates are manufactured to match individual patient CT data.

Surgeon-Specific plate designs include customizing plate thickness, screw hole placement, contour, and more.

AccuPlate® Optimized by AccuPlan® to work with complementary surgical planning and guides.

Save Time in the OR by eliminating the need to modify fixation plates during surgery.

Implant-Grade Titanium is highly resistant to corrosion and breaking with its excellent material performance and durability.

Milled Plates preserve strength and accuracy by eliminating mechanical stress induced by plate bending in the OR.

Process

START

Service request and **CT data** submitted via MedCAD.com/Start.

DESIGN

Surgeon AccuPlate® design is processed by MedCAD.

APPROVE

Surgeon approves AccuPlate® design for manufacturing or requests changes.

MAKE

Plates are CNC milled by MedCAD.

RECEIVE

Get AccuPlate® Patient-Specific Plate.

Required Information

CT Data

Service Request

Online Planning Session

To Order:
MedCAD.com/Start

Shipment:
MedCAD
501 S 2nd Ave,
Suite A-1000
Dallas, TX 75226

Digital File Transfer:
MedCAD.com/upload

Phone: +1 (214) 453-8864



This document is intended for physicians and healthcare professionals. Surgeons must always rely on their clinical judgement and use discretion when selecting products to use with an individual patient. MedCAD does not offer clinical advice and recommends that all surgeons are trained prior to use of any product in surgery. Always refer to the instructions for use, package inserts, and labels, including sterilization instructions, before using MedCAD products. © Copyright 2021, MedCAD, a DBA of VanDuzen, Inc. MedCAD and MC logos, AccuPlan, AccuPlate, AccuModel, and Custom Surgical Solutions are all trademarks of MedCAD.

Made in USA | PLT-PL-014, Rev A (193528)